care-e-vac

Medical Suction Equipment

Directions for Use and Maintenance

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Français

1. Cautions and Warnings

CAUTION: Federal law (USA) restricts this device to sale by or on the order of a physician (CFR 21 801.109).

- The care-e-vac unit must be used by trained personnel in the use of medical equipment.
- Read and follow the instructions of this guide.

Risk of explosion or fire: Do not use the machine near of flammable or explosive liquids or gases.

Risk of electric shock: Do not submerge this machine under any liquid. It could cause an electric discharge hurting to the people or damage the machine itself. Unplug the equipment before cleaning the unit. Use a minimum quantity of water or other liquids to avoid electric discharge.

Do not remove cover. Refer servicing to qualified personnel.

To avoid the risk of electric shock, this equipment must only be connected to a supply main with protective earth

Risk in extreme environmental conditions: Do not use this machine under environmental conditions that exceed the specified limits: -12 - 50 °C (10 - 122 °F) and 10 - 95% humidity; this could be dangerous to the security and affect the functioning in a negative way.

Risk of liquid entrance inside the pump: Do not start the machine if there is any deterioration in the overfilling system of the collector jar. Do not start the machine if there exists or if you suspect the existence of contamination owing to excess liquid. Return the device for service

Risk of contaminating the patient: Always use a clean collector jar and a new bacteria filter for each patient. After each use, it is necessary to disinfect the whole collector jar (see *Chapter 7.2*)

Risk for the machine due to a wrong connection: While changing the battery, make sure not to exchange the polarity, the poles have the same size and it is possible to connect them incorrectly (see *Chapter 7.1.2*).

- Remove the battery when the device is not used for long periods of time.
- To ensure proper grounding reliability, connect to “Hospital grade” receptacle or equivalent.
- Repairs are to be carried out by personnel qualified by Ohio Medical or one of its authorized dealers. This device must be clean of all blood or other organic materials prior to returning for service.
- Do not modify this equipment without authorization of the manufacturer.
- Please follow the instructions on the use of cleaning products and disinfectants.
- Please use only the accessories provided by Ohio Medical or its authorized distributors.
- Old electrical and electronic equipment must be disposed separately and may not be included in regular domestic waste. Disposal of device, components and accessories must strictly conform to local laws and regulations as set out by the relevant authorities.
- The configuration of the components of the product may be subject to variations without previous notification. We recommend that you consult with the Ohio Medical dealer to find out the latest version.

2. Introduction

2.1. Indications for use

The care-e-vac unit is an electrically powered portable suction equipment, intended for use in clinical institutions; the care-e-vac *PLUS* and care-e-vac *EMS* models are specifically designed for transport and emergency field use.

It is indicated for taking out secretions, blood or vomits from the respiratory track of the patients to let air flow freely to the lungs.

2.2. Contraindications

None known

2.3. Conformity with product standards and regulations

ISO 10079-1

ISO 10079-4

IEC 60601-1

IEC 60601-1-2

IEC 60601-1-6

IEC 60601-1-12

EN 1789

REGULATION 10 Rev.5 (E/ECE/324/Add.9/Rev.5) Uniform provisions concerning the approval of vehicles with regard to electromagnetic compatibility


3. Symbols

	Degree of protection provided by enclosures against dripping water		Distributor
	Class II (IEC 60601-1)		Date of manufacturing Year and month (yyyy-mm)
	Type BF (IEC 60601-1)		Reference number of the device
	Direct current (DC)		Serial number
	Attention. Consult accompanying documents		


4. Description of the equipment

4.1. General description

- A. Main unit (casing)
- B. Liquids collection jar
- C. Suction tubes (must be changed for each patient).
- D. Bacteria filter (must be changed for each patient).
- E. Electrical connection cables (depending of the models)
- F. ON/OFF switch
- G. Vacuum gauge
- H. Vacuum regulation control knob
- I. LED indicators

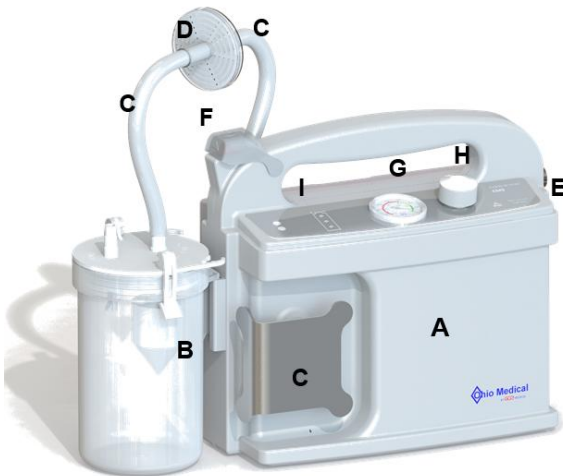4.2. Led indicators and vacuum control knob

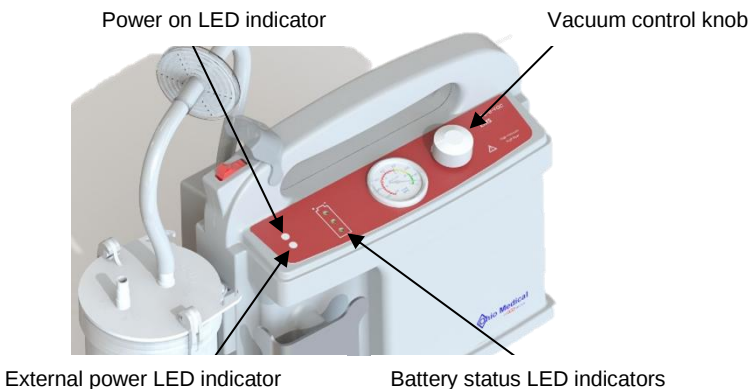Power ON LED indicator:

This green LED is continuously lit while the care-e-vac unit is switched ON.

External power LED indicator:

This green LED is continuously lit while external AC or DC-power is connected.

Battery status LED indicators:

3 LEDs: fully charged, 75-100% battery capacity

2 LEDs: 50-75% battery capacity


1 LED: 25-50% battery capacity

Device settings, state of battery, temperature, etc can influence battery status

Vacuum control knob: Turn right: increase vacuum
 Turn left: decrease vacuum

4.3. Electrical connection cables

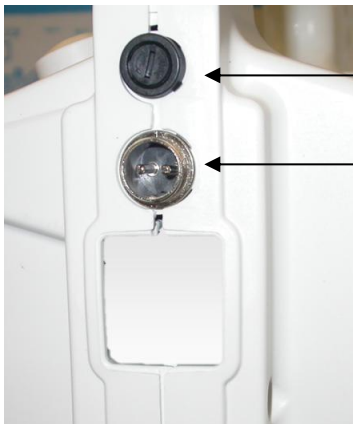12 VDC Fuse

12 VDC connector socket

5. Preparing of the equipment for its use

The care-e-vac suction unit is supplied assembled with the exception of the suction tubes to connect the bacteria filter between the suction unit and the collection jar; afterwards the equipment will be ready to be used. (see Chapter 5.1.).

The collector jar, the suction tubes and the bacteria filter are the only components that have to be changed for each patient and are easily removable. On assembling the equipment, special care has to be taken to respect the direction of the suctioning flow that is identified and recorded in the collector jar.

Never operate the unit without bacterial filter, the unit could be damaged

Do not forget to check the correct operation voltage

5.1. Previous assembly procedure

1. Connect the 5 cm (2") tube between the bacteria filter and the input of the suction unit.
2. Connect the 15 cm (6") tube between the bacteria filter and the outlet *VACUUM* of the collection jar.
3. Connect the suction tube (1.5mx8mm / 5ft x 0.3") between the outlet *PATIENT* of the collection jar and the suction terminal and roll it in the tube location of the unit.
4. Connect the electric power cable as required by the model.
5. Switch on the unit (ON/OFF switch).
6. Obstruct the flow by kinking the suction tube and checking the vacuum gauge to test that the unit makes vacuum without any leakage and that this level varies when the vacuum regulator is adjusted.
7. (Only in models with battery) Check the battery level while unit is running to charge battery, leave it connected to a power supply during the time needed (the full charge of a new battery that is run down is five hours).
8. Switch off the unit (ON/OFF switch).

5.2. Final assembling for a specific patient

1. Connect the suction device to a suction catheter specific for this application and patient.
2. Switch-on the unit (ON/OFF SWITCH).
3. Obstruct the flow in the suction catheter and adjust the level of vacuum prescribed in the vacuum gauge turning the vacuum regulator control knob.
4. The unit is ready to put the suction catheter to the patient.

6. Functioning and Operation

-  Do not start the machine if you suspect the presence of explosive or inflammable liquids or gases, water inside the machine or environmental conditions that are not within the specified limits.
-  Do not start the unit if the liquid level of the jar exceeds 50% of its capacity.
-  When operating unit in battery mode, the user must monitor the LED's to ensure enough battery power is left to finish procedure.
-  MR unsafe

- Verify the assembly and test it as detailed above in *Chapter 5*.
- Proceed to a visual inspection: check that there are no one damaged components and the unit and the collection jar are clean.
- Place the unit on a flat surface.
- Apply the suction catheter to the patient.
- When the liquid level reaches 50% of the total volume of the collection jar, it is recommended to empty or change it before the safety valve acts, in this way the risk of the entrance of humidity in the bacteria filter is minimized.
- Once used, disassemble suction line and proceed to its cleaning / substitution as indicated in *Chapter 7*. After this will be necessary to assemble the unit as detailed in *Chapter 5*.

7. Maintenance and Cleaning

7.1. Periodic (planned) maintenance

Clean the unit after every use or maintenance.

Perform an operation test before each use and after any maintenance operation and at least once a month (see *Chapter 5.2*).

Perform a battery capacity test every 6 months (see *Chapter 7.1.1*) and charge the battery for 5 hours after each use connecting it to the mains.

It is recommended to replace the battery every 3 years, or before in the case of any failure in the scheduled tests (see *Chapter 7.1.2*).

7.1.1. Test of capacity of the battery

1. Charge the battery for 6 hours and disconnect from the supply.
2. Switch-on the unit and turn the vacuum regulator knob to high vacuum, make sure that the vacuum gauge is at the correct level.
3. Suction speed: Measure the time needed to suck 500 ml of water.
4. Repeat this step 5 times; the average time should be less than 10 seconds.
5. Let the unit work continuously with free flow for 20 minutes.
6. If the unit stops before 20 minutes or it does not pass the suction speed test, substitute the battery with a new one pursuant Chapter 7.1.2.

7.1.2. Battery replacement

1. Unplug the unit from the mains.
2. Remove the 6 screws from the back side.
3. Maintain the unit in vertical position and separate the both parts of the casing.
4. Remove the battery and recycle it pursuant the local rules.
5. Install the new battery with its terminals in its correct polarity.
6. Close the casing and install the 6 screws in the back-side.
7. Charge the battery for 5 hours and test its capacity (*Chapter 7.1.1*).


7.1.3. Fuse replacement

Unplug the unit from the mains

Replace the fuse with a new one of the same rating. (See Technical specifications)

12/24Vdc fuse:
turn left using a
flat head
screwdriver and
take the fuse
holder out

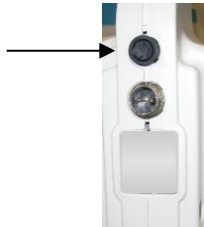7.2. Cleaning and Disinfection

Disconnect the machine from any external mains before cleaning it. Use a minimum quantity of water or other liquids to avoid an electric discharge.

Contaminated devices must be thoroughly cleaned prior to disinfection or sterilization

Personal protective equipment must be worn when handling contaminated equipment, particularly tubing, filter and collection jar (protective clothing, disposable latex gloves, mask, eye protection)

If any liquid accidentally enters the main unit or the pump, the unit may be damaged. For this reason, we recommend not to use the aspirator itself as a method to clean the collection jar (taking in water with disinfectants) nor connecting the suctioning tubes in the opposite way (the safety device would not work). The failure should avoid the warranty and the technical specifications under the manufacturer's responsibility.

In the event that contaminated fluids have moved past the internal unit, remove the tubing and the filter, clean the casing (7.2.1), and decontaminate the collection jars (7.2.2). Return the device for service.

Reusable polycarbonate jars can be reused up to 40 times and Polysulfone jars up to 100 times. Before each use, check that collection jars are not damaged.

To obtain additional information including questions on infection control procedures, please contact Ohio Medical at ohiomedical.com or 866-549-6446.

7.2.1. Cleaning of the main unit (casing)

1. Use a cloth or sponge that is dampened with a mild detergent to clean the external surfaces of the main unit.
2. Repeat the previous step with a mild disinfectant compatible with ABS and PC.
3. Use a cloth or sponge dampened with water and wipe the surfaces again. Dry the surfaces using a clean cloth or a paper towel.

 Never use dissolvent or abrasive products.

7.2.2. Decontamination of the liquids collection jar

-  It is the user's responsibility to qualify any deviations from the recommended method of processing.
-  Follow the decontamination product's instructions. Check the decontamination product compatibility with the materials
-  After decontamination, label products with the date
-  The disposable jars should be disposed according to the local rules.


After each use, decontaminate the collector jar pursuant the next procedure:

A. Disconnect the 3 tubes and the filter and replace them for new ones.

B. Disassemble the cover and the floating device from the transparent jar:

Free the two locking clamps to open the jar, lift the cover and pull down the floating device.

If dirt is observed in the rubber quad-rings of the floating device or in the locking clamps, they can be disassembled (it is recommended to warm up in water at 50°C (120 °F) to expand the rubber rings).

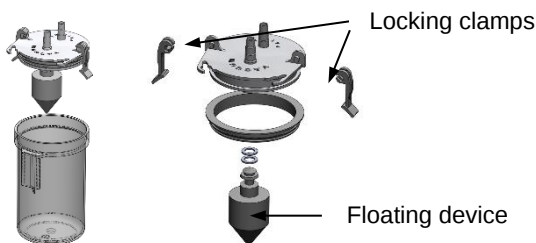C. Washing and rinsing (manual or mechanical):

Wet in lukewarm water (<43°C / <110 °F) before adding any detergent; Wash parts thoroughly in warm water using a compatible detergent (rub with brushes, if necessary); Rinse thoroughly with sterile water, by immersing it completely in three separate copious volumes of water. Each rinse should be a minimum of one minute. Do not reuse the water. Finally, allow to dry.

- Alkaline or neutral detergents: non-ionic tension-actives concentration lower than 15%, or ionic tension-actives concentration lower than 30%.
- Enzymatic detergents: enzymes concentration lower than 15%.

D. Disinfection / Sterilization:

Disinfected/Sterilized jars should be dated and stored in clean and dry areas to avoid recontamination

Avoid direct contact with the autoclave chamber, because the material can be damaged. It is recommended to use a basket.

All parts can be processed pursuant one of the next 2 methods:

- **Steam Sterilization in gravity autoclave:**

Polycarbonate (PC) jars: (max. 40 times) at 121°C (250 °F) and 1.3 kg/cm² (18,5 psi) during 25 minutes, drying time: 30 minutes.

Polysulfone (PSU) jars: (max.100 times) at 132-136°C (270-277 °F) and 2.4 kg/cm² (34 psi) 15 minutes, drying time: 30 minutes.

- **High level Disinfection:**

With Glutaraldehyde (i.e. *Cidex*): 2% solution at ambient temperature during 60 minutes. Rinse thoroughly with sterile water, by immersing it completely in three separate copious volumes of water. Each rinse should be a minimum of one minute. Do not reuse the water. Finally, allow to dry.

E. Inspection, reassembly and perform the function test pursuant Chapter 5.

8. Troubleshooting information

FAULT	CONDITION	ACTION
The unit does not operate	Control Unit not turned on	Turn Main Switch "ON/OFF"
	Wrong operating voltage	Check mains voltage output
	Mains not connected	Connect mains to the Unit
	External power indicator lit	Send the unit to repair
	Defective Fuse	Check fuse and replace it
The unit does not work from its battery	Battery status LEDs are not lit	Place the unit on charge. If does not operate after charging, replace the battery
The unit operates, but with little flow or little vacuum	Floater valve is closed	Unplug the Filter/Vacuum tubing to release the vacuum
	Collection jar full	Remove and replace the jar
	Defective connection from the Unit to the jar	Install the Vacuum Connector Tubing correctly
	Patient Suction Tubing twisted or blocked	Untwist, clear or replace the tubing

9. EMC Information

Care-e-vac equipment meets the IEC 60601-1-2 EMC requirements for medical devices.

Medical electrical equipment needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information. Portable and mobile RF communications equipment can affect medical electrical equipment.

10. Technical Specifications

	care-e-vac <i>PLUS</i>	care-e-vac <i>EMS</i>
RANGE	High Vacuum High Flow	
SUPPLY	12 Vdc Battery	
VACUUM	84 kPa 630 mmHg	
FLOW	30 L/min	
POWER	85 W	
BATTERY	12 V 4,5 Ah. Ac-Pb	
AUTONOMY (battery)	65 min	
BATTERY CHARGER	YES	External (optional accessory)
10 g tested wall bracket	-	YES
WEIGHT	4,2 kg (9,3 lbs)	
Fuse	DC F8AL250V	
Dimensions (packaged)	460 x 210 x 290 mm (18 x 8,3 x 11,4 in)	
Noise level	60 ± 1,5dB (AS) / 1 m	
Storage Conditions	-4°F to 140°F (-20°C to 60°C) 10% to 95% RH 70 to 106 kPa	
Operation Conditions	-10°F to 122°F (-12°C to 50°C) 10 % to 95 % RH 70 to 106 kPa	

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