





## **Extremely resilient and safe**

Our innovative HERZmed PROFI transfer sheet is a waterproof single-use sheet for transferring patients for use in rescue services as well as in hospitals and nursing homes. The load-bearing and absorbent capacity of the transfer sheet has been tested by DEKRA (DEKRA Report, Report No. 340785800).

The sheet is extremely robust, even after absorbing liquid. The integrated, high-performance absorption zone (three layers consisting of PE films, cellulose fluff and fleece) prevents unwanted waterlogging - to the benefit of the patient and the facility.


#### **Details**







Suction area Backside Single packaging

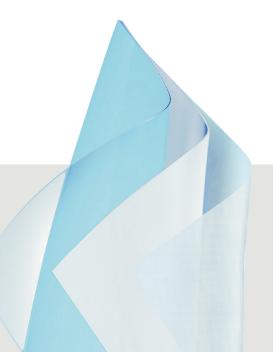


## **Highlights**

#### THREE-LAYER SYSTEM:

- 1. PE-film: This material prevents liquids from leaking through. The roughened blue side avoids slipping of the cloth.
- 2. Cellulose fluff: Due to the special rhombic-shaped arrangement, the escaping liquid is quickly distributed in the absorption body.
- 3. Fleece: The upper side consists of a skin-friendly fleece fabric. The fleece layer forms a comfortable patient support with a pleasant feel.



#### **Extremely resilient**

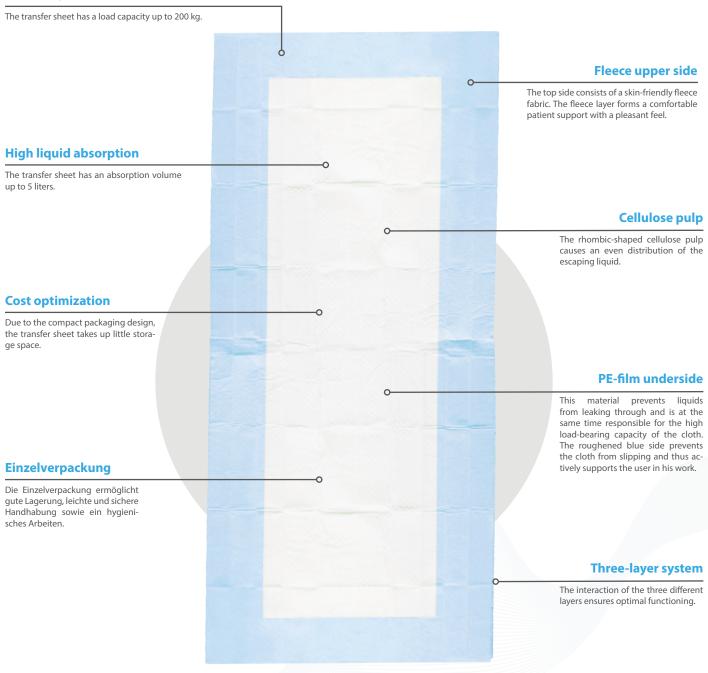

# Hygienic working environment

Due to the generous dimensions, the entire work area is safely covered and effectively protected against damaging influences caused by moisture. The transfer sheet Profi can minimize so-called cross-infections for patients and staff.

This is a single-use product!



## **Technical specifications**

Size:

Main fabric: 220 x 110 cm (L x B) Suction zone: 200 x 50 cm (L x B)

**Weight:** Ca. 350 g

Coloring:

**Packaging:** Single packed

**Absorption volume:**Absorption: Up to 5 liters

**Load capacity:** Up to 200 kg (even after liquid absorption) Material:

Under layer: Polyethylene (PE) Middle layer: Cellulose fluff

Sterility:

### **Order information**

ArtNo.	Description
V100048	Transfer sheet profi, Model 2020

Order and advice by:				

**HERZ**med GmbH

Herkulan-Schwaiger-Gasse 16 82487 Oberammergau

Tel.: +49 (0) 911/14 88 959 - 0 Fax: +49 (0) 911/14 88 959 - 9 Mail: herzmed@herzmed.de

Web: www.herzmed.de

The products shown are only an excerpt of our product range. Exact details and product information are available on request or at <a href="https://www.herzmed.de">www.herzmed.de</a>. All prices marked with an asterisk are exclusive of VAT. All brand names used remain the property of the respective brand owner. We assume no liability for errors, possible misprints and color deviations. Our general terms and conditions apply, which you can also view at <a href="https://www.herzmed.de">www.herzmed.de</a>.

